

Installation Instructions



9A302001 Welch Allyn Light Bracket (Use w/ 48700 Light)

Intended Use

To mount a Welch Allyn 48700 Light to table.

Restrictions

Requires Base Rail (see accy. # below):

9A184 [Models 311, 319, & 75(-015)]


9A205 [Models 411(-016 thru -018), & 75(-016 / -017)]

9A272 [Model 419(-003 thru -006)]

Kit Components


<u>Description</u>	<u>Part Number</u>	<u>Qty</u>
Bracket	9A302001	1

Installation (1)

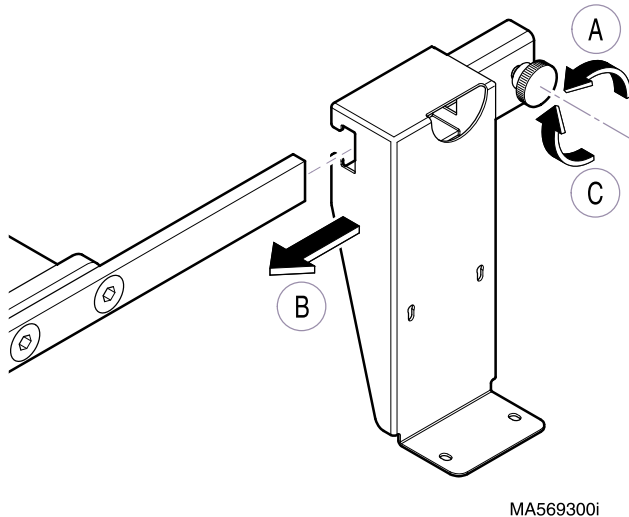

- A. Loosen knob.
- B. Slide hanger assembly onto base rail.
- C. Tighten knob.

Installation (2)

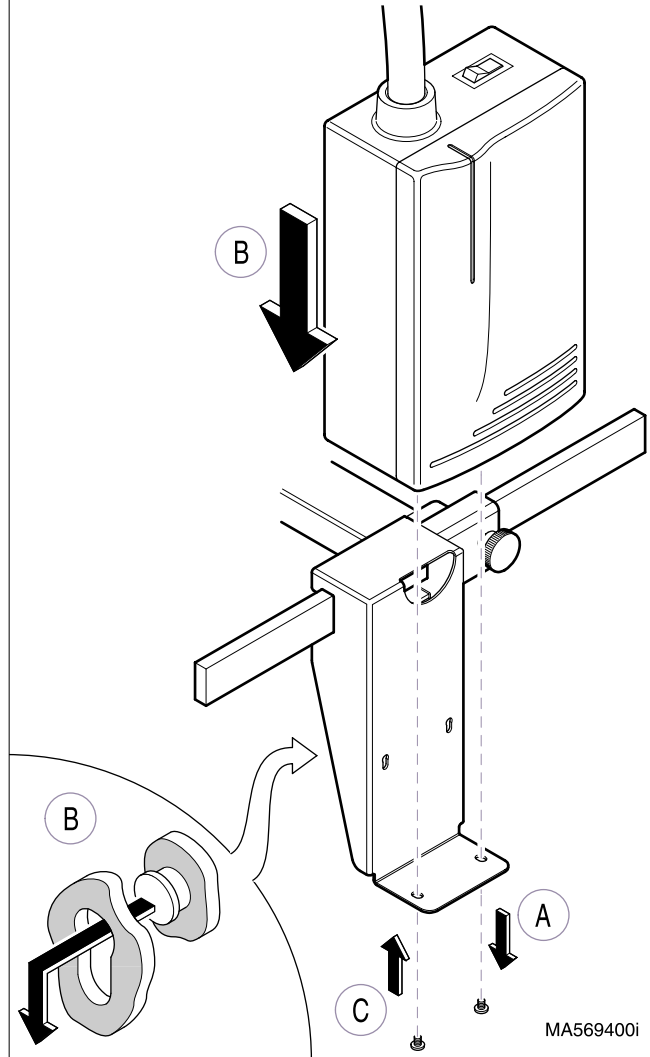

- A. Remove two screws from bottom of light.
- B. Slide light onto bracket until pins lock into slots in bracket.
- C. Secure light to bracket with two screws removed in step A.